

SUMMARY

shipped at RT; store at 4 °C

For research use only

Product

pUC18 high copy cloning vector for replication in *E. coli*, suitable for “blue-white screening” technique.


Description

pUC18 is a small, high copy cloning vector for replication in *E. coli*. It has been constructed using the ampicillin resistance gene and the pMB1 origin of replication from pBR322. The pMB1 of pUC18 differs from the pBR322 origin by a single point mutation and the lack of the *rop* gene, leading to a high copy number. Additionally, pUC18 contains the *lac* operon of *E. coli* with CAP binding site, *lac* promoter (P_{lac}), Lac repressor (LacR) binding site, and the 5'-terminal part of the *lacZ* gene encoding for the N-terminal part of β -galactosidase (source – M13mp18 phage vector). This 5'-terminal part of the *lacZ* gene contains the multiple cloning site (MCS), and its expression is IPTG inducible. It is capable of intra-allelic α -complementation of a partial deleted chromosomal *lacZ* copy (*E. coli* host strain: *lacZ* Δ M15, e.g., DH5 α , DH10B, JM101, JM109). In the presence of IPTG, transformants expressing both fragments of the β -galactosidase (the vector encoded N-terminal part and the chromosomal encoded C-terminal part) will form a functional enzyme and can be detected as blue colonies on agar plates containing X-Gal. Cloning into the multiple cloning site will lead to a nonfunctional N-terminal fragment of the β -galactosidase and to the abolishment of α -complementation. White colonies will form on X-Gal/IPTG plates.

Cloning Vector pUC18

Product Information Sheet
V33002

Vector Map

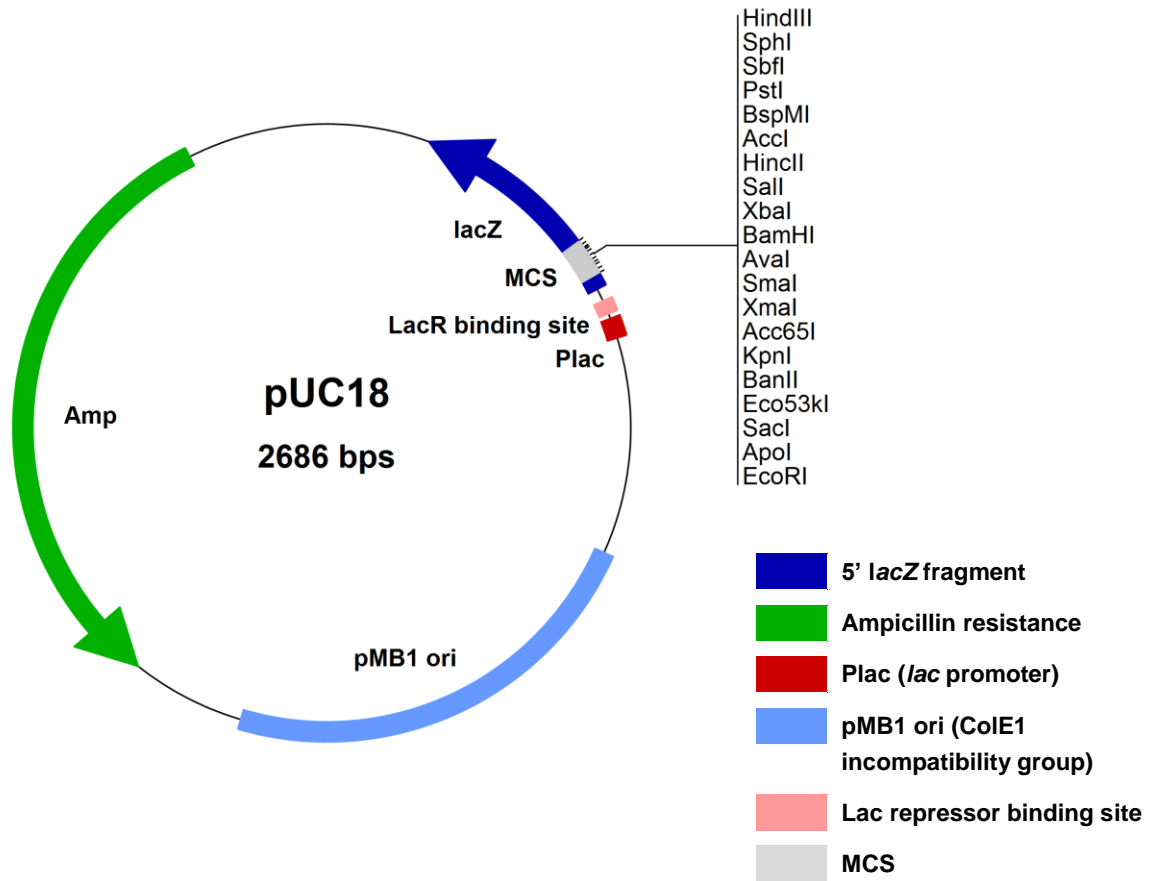

Quality Warranty

DNA concentration and purity was checked by UV spectrophotometry. All restriction sites specified in the vector map were checked by sequencing. Functionality of α -complementation was checked by transformation and plating the transformants on IPTG/X-Gal agar plates.

Cloning Vector pUC18

Product Information Sheet
V33002



References

Yanisch-Perron C (1985) Improved M13 phage cloning vectors and host strains: nucleotide sequences of the M13mp18 and pUC19 vectors; Gene 33, 103-119

Order Information, Shipping and Storage

Order#	Product	Amount
V33002	pUC18, lyophilized DNA	25 µg
shipped at room temperature (RT); store at 4 °C. Once the DNA has been dissolved in sterile water or buffer we recommend storage at -20 °C.		

Related Products

Order#	Product	Amount
MTAQK0	MoBiTaq-K (25 U/µl)	250 U
STAQ02	SuperTaq (5 U/µl)	250 U
STAQH1	Super Taq-HC (15 U/µl)	250 U
ENZ-286-1PS	Recombinant T4 DNA Ligase	20,000 U
GE-TLK0110-1	TurboLigation™ Kit	100 rxn
V33202	pUC19 vector DNA	25 µg
RIBA25	RNAse A, 90 U/mg (Kunitz)	25 mg
A1414-25GMAG	Ampicillin, sodium salt	25 g
I1312-1gAG	IPTG (Isopropyl-Beta-D-thiogalactoside)	1 g
X1015-5gAG	X-GAL	5 g
04004G	MoBiTec Agarose LE	500 g

Contact and Support

MoBiTec GmbH ● Lotzestrasse 22a ● D-37083 Goettingen ● Germany

Customer Service – General inquiries & orders

phone: +49 (0)551 707 22 0
fax: +49 (0)551 707 22 22
e-mail: order@mobitec.com

Technical Service – Product information

phone: +49 (0)551 707 22 70
fax: +49 (0)551 707 22 77
e-mail: info@mobitec.com

MoBiTec in your area: Find your local distributor at www.mobitec.com