


Anti-SPINK1 (hPSTI) antibody, monoclonal mouse

Product Information Sheet
PSKAN2

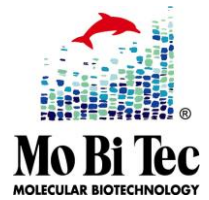

SUMMARY

shipped on blue ice; store at -20 °C

For research use only

monoclonal mouse IgG, (Clone R XXIII)

Introduction

Antibody to human pancreatic trypsin inhibitor domain (hPSTI); recognizes soluble hPSTI as well as hPSTI in the pSKAN phagemid pIII-hPSTI fusion protein.

Purification: Protein A purified

Working Dilution: 1:250 to 1:1000 for immuno-blotting (Western Blots) and ELISA; for incubation the blots should be welded in plastic foil, reducing the assay volume to a few ml.

Amount: 100 µl (1.66 mg/ml)

Storage Buffer: PBS, 0.9 mg/ml NaN₃; pH 7.4

References

C. Férec *et al.* (2007) *Eur J Hum Genet.* 15(9):936-42. PMID 17568390
M.Tesar *et al.* (1995) *Immunotechnology* 1,53-64. PMID 9672201

Order Information, Shipping and Storage

Order#	Product	Quantity
PSKAN2	Anti-SPINK1 (hPSTI) antibody, monoclonal mouse	100 µl
shipped on blue ice; store at -20 °C		

Contact and Support

MoBiTec GmbH ● Lotzestrasse 22a ● D-37083 Goettingen ● Germany

Customer Service – General inquiries & orders information

phone: +49 (0)551 707 22 0
fax: +49 (0)551 707 22 22
e-mail: order@mobitec.com

Technical Service – Product

phone: +49 (0)551 707 22 70
fax: +49 (0)551 707 22 77
e-mail: info@mobitec.com

MoBiTec in your area: Find your local distributor at www.mobitec.com