

**ProFoldin**

10 Technology Drive, Suite 40, Number 188

Hudson, MA 01749-2791 USA

Phone: (508) 735-2539

FAX: (508) 845-9258

www.profoldin.cominfo@profoldin.com

INSTRUCTIONS

ProFoldin *S. pneumoniae* DNA Helicase

S. pneumoniae* DNA Helicase – for 100 assays*Catalog No. DNAB-100SP**

Protein construct: Wild-type *S. pneumoniae* DNA helicase purified from a bacterial expression system.


MW: 50 kDa

Enzyme concentration: 20 μ M

Enzyme activity assay: The ATPase activity of DNA helicase is measured by using the **DNA Helicase ATPase assay Kit** (Catalog No. DNAB100K).

Storage temperature: -20 or -80°C. Do not freeze-and-thaw repeatedly.

Enzyme dilution: Use the 1 x assay to dilute the enzyme just before the assay. Do not store diluted enzyme solution

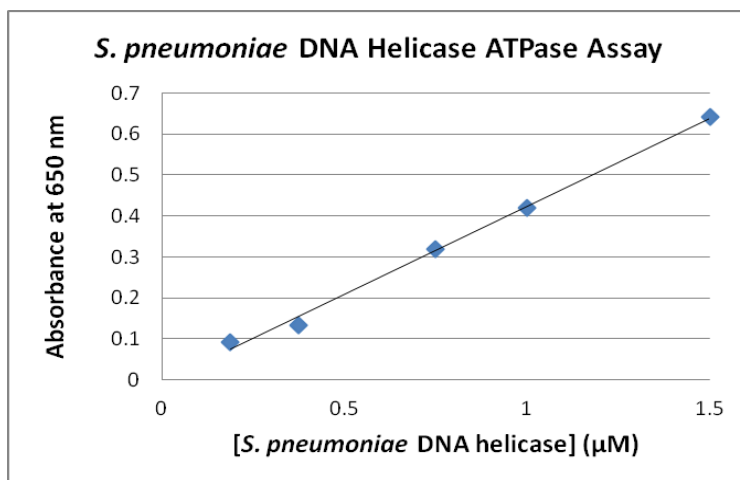

The ***S. pneumoniae* DNA Helicase – for 100 assays** (Catalog No. DNAB-100SP) includes 33 μ l of 100 x *S. pneumoniae* DNA helicase (20 μ M). It is for 100 assays.

Assay Protocol using the DNA Helicase ATPase assay Kit

1. Reagent preparation:

For each 10 assay reactions,

- (1) Prepare 297 μ l of premix composed of 261 μ l of H₂O, 33 μ l of 10 x Buffer and 3.3 μ l of 100 x *S. pneumoniae* DNA helicase.
-



ProFoldin

10 Technology Drive, Suite 40, Number 188

Hudson, MA 01749-2791 USA

Phone: (508) 735-2539

FAX: (508) 845-9258

www.profoldin.com

info@profoldin.com

INSTRUCTIONS

(2) Prepare 33 μ l of 10 x Enzyme substrate by mixing 3.3 μ l of 100 x ATP and 3.3 μ l of 100 x DNA and 26.4 μ l of water.

2. Reaction:

Mix 27 μ l of the premix with 3 μ l of the 10 x Enzyme substrate in each well. Incubate the reaction mixture at 37°C for 60 min.

3. Detection:

Add 45 μ l of the Dye MPA3000 into the 30 μ l of the reaction mixture. Incubate for 5 min. Measure the light absorbance at 650 nm.

Note: The final concentrations for the ATPase assays of the helicases are 20 mM HEPES, pH 7.5, 20 mM potassium glutamate, 1 mM DTT, 0.005% Triton X-100, 10 mM MgCl₂, 20 μ g /ml DNA, 0.25 mM ATP and 200 nM DNA helicase. A negative control reaction can be the reaction mixture without addition of ATP or enzyme.

Assay optimization for enzyme inhibition

The assay can be optimized in terms of assay window, assay linearity and sensitivity to competitive inhibitors. ProFoldin offers HTS assay development service. For more information, please visit our website at <http://www.profoldin.com/services.html>.

Reference

Nakano T., et al, Translocation and Stability of Replicative DNA Helicases upon Encountering DNA-Protein Cross-links, J. Biol. Chem. 288: 4649-4658 (2013).

Related products

DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100K

E. coli DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100KE

H. influenzae DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100KH

S. pneumoniae DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100KN

P. aeruginosa DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100KP

S. aureus DNA Helicase ATPase assay Kit Plus-100

Catalog No. DNAB100KS

E. coli DNA Primase Assay Kit Plus

Catalog No. EGA100KE

H. influenzae DNA Primase Assay Kit Plus

Catalog No. HGA100KE

S. pneumoniae DNA Primase Assay Kit Plus

Catalog No. PGA100KE

S. aureus DNA Primase Assay Kit Plus

Catalog No. AGA100KE

For more information of drug targets and enzyme assays, please visit www.profoldin.com or send emails to info@profoldin.com.



ProFoldin

10 Technology Drive, Suite 40, Number 188

Hudson, MA 01749-2791 USA

Phone: (508) 735-2539

FAX: (508) 845-9258

www.profoldin.com

info@profoldin.com

INSTRUCTIONS

ProFoldin. Revised 07/2017. All rights reserved.