


MagWell™ 96-Well Magnetic Separator

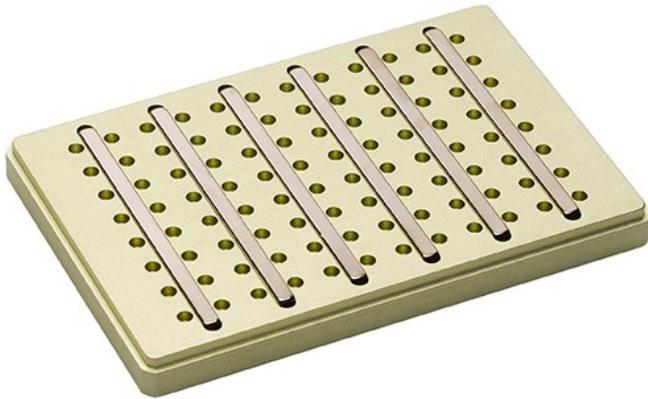

MagWell™ 96-Well Magnetic Separator features a unique recessed well design with bar magnets imbedded in an aluminum block that provides significant performance advantages in magnetic separation for 96-well applications.

Benefits

- **Efficient Particle Washing** - vertical orientation of the bar magnetic field attracts and suspends the magnetic particles both to the side and up from the bottom of the wells. Pipet tips may be inserted to the bottom of the well, insuring thorough and efficient washing and removal of fluids
- **Minimal Particle Loss** - suspended particles held up and to the side from the point of fluid aspiration minimizes particle loss
- **Maximum Sample Concentration** - samples can be efficiently recovered in 5 μ l vs the recommended 40 μ l for a typical donut style separator
- **Universal Compatibility** - compatible with all standard microtiter plates